

Bariatric patients



Liko UltraSling, Mod. 06

With a 500 kg. capacity, UltraSling is designed for the heaviest bariatric patients. UltraSling is used together with Liko's UltraTwin lifting system and is available in different sizes. With the UltraTwin lifting system, the patient can easily be positioned while in the sling during lift and transfer.



ComfortSling Plus, Mod. 300, 310

ComfortSling Plus is an upgraded version of ComfortSling (Mod. 30, 31) with improved comfort and fit. ComfortSling Plus enables an upright seated posture and provides support for the entire back. ComfortSling Plus is applied and removed while the patient is lying down, and is then left in the chair after lifting. Model 310 features a hygiene opening.

Basic models for safe and easy lifting



Mod 10

Liko OriginalSling, Mod. 10, (11)

OriginalSling has a semi-high back that allows the patient to keep one or both arms outside the sling.

OriginalSling is appropriate for the active patient, and it is easy to work with.



Mod 00

Liko UniversalSling, Mod. 00, 02

UniversalSling is a basic model that is suitable for most patients. UniversalSling is ideal for elderly patients with limited mobility.

It allows an upright sitting posture with support for the entire back.

Mod. 02 enables extra lifting height.



Mod 20/26

OriginalSling /Soft Original Sling HighBack, Mod. 20, 26

OriginalSling/Soft OriginalSling HighBack enables a semi-reclining seated posture, with support for the head and is therefore ideal for patients with poor torso stability. Soft OriginalSling HighBack has been specially adapted for improved fit and comfort, with comfortable head support in the seated position. OriginalSling/Soft OriginalSling is very suitable for lifting to/from the floor.



Mod 350/360

ComfortSling Plus HighBack, Mod. 350, 360

ComfortSling Plus HighBack is an upgraded version of ComfortSling (Mod. 35, 36) with improved comfort and fit. ComfortSling Plus HighBack enables a semi-reclining seated posture and provides support for the entire back. It features a soft head support. ComfortSling Plus HighBack is applied and removed while the patient is lying down, and is then left in the chair after lifting. Model 360 features a hygiene opening.



Mod 00

UniversalSling is perfect for lifting with leg supports overlapped under both legs.



Mod 22

Liko SilhouetteSling, Mod. 22

SilhouetteSling enables a semi-reclining sitting posture with comfortable head support. It has been specially developed for transfers to close-fitting wheelchairs where the sling can be left after transfer to the chair is completed. SilhouetteSling is thin and pliable, without handles or bulky parts.

Lifting to the toilet

Models that enable dressing/undressing during the actual lifting operation.



Mod 41



Mod 45/46

HygieneSling, Mod. 40, 41, 45, 46

The HygieneSling is an easy-to-use sling that facilitates dressing and undressing during lifting to and from the toilet. The HygieneSling is appropriate for patients with a certain degree of torso stability.

HygieneSling Mod. 41 is equipped with a belt and has improved fit.

HygieneSling Mod. 45 is equipped with a SafetyBelt and lifts around the upper body for extra security.

Model 46 is equipped with a SafetyBelt and an extra back support.



Mod 50



Mod 55

Liko HygieneVest, Mod. 50, 55

HygieneVest is a unique sling that is specially designed to be safe and secure for lifting patients with poor muscle tone to and from the toilet.

HygieneVest enables an upright sitting posture with firm support around the upper body while leaving much of the lower body free for dressing and undressing.

HygieneVest HighBack Mod. 55 is equipped with a head support.

Disposable Slings Solo



Mod 25

Solo HighBack Sling, Mod 25

To meet the need for slings for individual patients and as a way to avoid laundering, Liko has developed a series of disposable slings. Solo HighBack is appropriate for most patients and can be used in all the most common lifting situations.

Amputees



Mod 70

Liko AmputeeSling, Mod. 70, 75

Liko AmputeeSling is specially created to provide full security, even when lifting people with high double-sided amputations.

The extra strap guide keeps the leg supports together for safety reasons.

Liko Amputee HighBack Sling Mod. 75 is equipped with a head support.

Lifts for support and safety in the standing position



Mod 60, 64



Mod 60



Mod 64

Liko MasterVest, Mod. 60 and 64

MasterVest makes it easier to lift a person who can support their weight in a standing position. MasterVest is an excellent aid to standing-, gait- and balance training. The lower body is free for dressing and undressing during transfers to and from the toilet.

Liko MasterVest, Mod. 64, with front fasteners, can be applied and removed by the patient. Leg harnesses are available as an accessory.



Mod 92



Mod 920

Liko LiftPants, Mod. 92, 920

Liko LiftPants facilitate safe, secure standing- and gait training, giving patients the confidence to take new steps. LiftPants enable freedom of movement, while relieving some of the burden of body weight. They lift safely, allowing the patient to move on his/her own without the risk of falling.

Sensitivity to pressure



Mod 10

For patients who are particularly sensitive to pressure, several padded models with synthetic sheepskin are available. The padding also has a reinforcement insert that can be removed if necessary.

Net polyester – an airy alternative



Mod 00

Net polyester fabric is a suitable material for a sling that can be left on the patient after lifting. The material is airy, yet strong, and several versions of slings in polyester net are available. Slings in net polyester are ideal for bathing and showering.

Bath and shower



Mod 00

For lifts to the bath and shower, there are several models in durable, non-absorbent plastic-coated net fabric that feels comfortable next to the skin, even when wet. Slings in plastic-coated net are ideal for bathing and showering.

Liko sling overview, adult sizes

[illegible]

- = Standard
- = Made-to-order

Lifting accessories



Repositioning Sheet

RepoSheet. Original is an excellent aid when a patient has to be moved higher up or turned in bed. The RepoSheet is bedded under the patient and can be used in combination with, or instead of, an ordinary bed sheet.



Sling Cross-bar

Sling Cross-bar enables a slightly more reclined sitting posture than a two-point slingbar and is ideal for use with many slings in Liko's product range. Two models are available.



Extension loops

Liko extension loops are intended for quick and easy adaptation of lifting accessories. There are two standard lengths: for 12 cm and 22 cm (4.7 or 8.7 inch.) extensions of the normal strap loop lengths. Available in green and grey.



Padded Cases


Padded Cases are available as an extra accessory for all our slings. They are fitted with synthetic sheepskin on one side and the sling's material on the other. Indicate which sling (model, size and material) the Padded Case is intended for when ordering.



LikoScale

A digital scale that can be used in combination with all of Liko's lifts. Connected above the slingbar, LikoScale is easy to use and easy to read. There are LikoScales for weighing patients up to 350 kg (770 lbs.).

Universal SlingBars

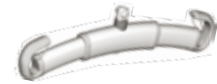

Universal SlingBar 350

This model is especially popular where lifting children is concerned.




Universal SlingBar 450

An allround model that is ideal in most cases.



Universal TwinBar 670

A slingbar for larger patients or in situations where greater width across the shoulders is desired.

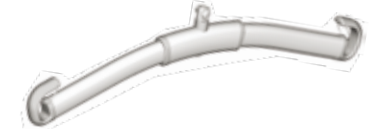

Universal SlingBar 600


This slingbar is mainly intended for use with ComfortSlings or when a sling is applied with separate leg supports "overlapped" under the legs.



Universal SideBars

Used to achieve a four-point suspension to enable a semi-sitting posture. Normally used in combination with Universal SlingBar 450 or 600.

See more accessories at www.liko.com

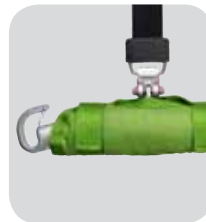

Hill-Rom

Enhancing Outcomes for Patients and Their Caregivers.™

France +33 (0)2 97 50 92 12
 United Kingdom +44 (0)1530 411000
 Deutschland +49 (0)211 16450-0
 Nederland +31 (0)347 32 35 32
 Italia +39 (0)2 950541
 Suisse/Schweiz +41 (0)21 706 21 30
 (deutschsprachig) +41 (0)21 706 21 38
 Österreich +43 (0)2243 28550

Ireland +353 (0)1 413 6005
 Spain +34 936 856 000
 Portugal +351 210 991 891
 Nordic Countries +46 (0)920 47 47 11

Export Contact your local Hill-Rom distributor or contact your Area Manager via website or call +1 812 934 8173

www.hill-rom.com
www.liko.com

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.
 © 2010 Hill-Rom SARL, ALL RIGHTS RESERVED

Lifting accessories.



Adult Slings

Hill-Rom

Enhancing Outcomes for Patients and Their Caregivers.™

Liko lifting accessories

It started with a sling made entirely of textile fabric. At the time, this was revolutionary. Now Liko has the world's largest assortment of slings and lifting accessories.

Together with our many varieties of slingbars in different widths and with different functions, slings can be used to lift patients in all types of situations.

For more than 30 years, we have been developing functional, easy-to-use lifting accessories in close collaboration with caregivers and patients. With growing demands for lifting and as knowledge in this area increases, Liko's development process continues.

Many models

All models have unique features to suit different lifting needs. Some are referred to as basic models, since they cover a wider range of lifting needs for different people. Other models are more specialized for optimal lifting in narrower application areas.

Children Slings

Children are not small adults. We have a large assortment of slings especially developed for children. Please refer to our website or our special brochure for children's accessories for more information.

Individual fitting and testing

To ensure maximum safety and comfort for patients and caregivers, each lifting accessory should be tested individually.

Sizes

The right choice of size is essential for safety, function and comfort.

Each model has individual sizes, as indicated in the table "Liko sling overview".

Instructions

A detailed instruction guide accompanies each sling. Everyone who uses the products **must** read and fully understand the instructions.

The instructions include recommendations concerning the most appropriate choice of slingbar for each sling.

Marking

Marking labels on the slings indicate the model, size and version as well as laundering instructions and max. load.

Inspection

Slings have a limited working life. Use and laundering eventually cause fabric and seams to wear. All lifting accessories must be inspected regularly. Never use damaged lifting accessories!

Liko's slings are tested by accredited test institutes and comply with all requirements for MDD Class I products and all requirements according to the harmonized EN ISO 10535 standard.